

CONTINUOUS DRYER TYPE «KTF»

THE FLEXIBLE CONTINUOUS HOT AIR DRYER

EFFICIENT APPLICATION

Continuous drying under production conditions with the «KTF» - everywhere, where the required sample length makes the efficient application of a large machine impossible.



POSSIBLE APPLICATIONS

The modular design of the «KTF» permits a variety of possible applications. Configured especially to your requirements, the «KTF» is used for the following working processes:

- Drying, curing and condensation in controlled hot air.
- Thermolising
- Drying and thermolising with two drying chambers in one double «KTF».

Continuous as well as discontinuous sequence of the working process with only one «KTF»:

- Continuous process flow in combination with MATHIS Padder type «HVF» or from roll to roll.
- Discontinuous process flow with specific material holders.

Expandable functionality of the «KTF» through complementing with process-focussed add-on components.

SCOPES OF APPLICATION

The scopes of application of the «KTF» are the textiles industry, the pharmaceutical industry, foil and paper manufacturers, but also research institutes and development departments of most branches.


PERFECTED FUNCTIONING PRINCIPLE WITH GREAT RESULTS

Continuous process flow in combination with MATHIS padder horizontal, vertical type «HVF»


The material sample is unrolled from the «HVF», treated process-related and guided through the «KTF» for drying. At the outlet of the «KTF» the dried material sample is rolled up automatically and is thus available for the next processing. The process speed is preset by the «HVF» in this configuration.

From roll to roll

The material sample to be dried is unrolled from the unwinding device of the «KTF» and guided through the dryer. At the outlet of the «KTF» the dried material sample is rolled up by an upwinding device, which is designed as a rising roll batcher. Optionally an upwinding with core or sleeve is available.

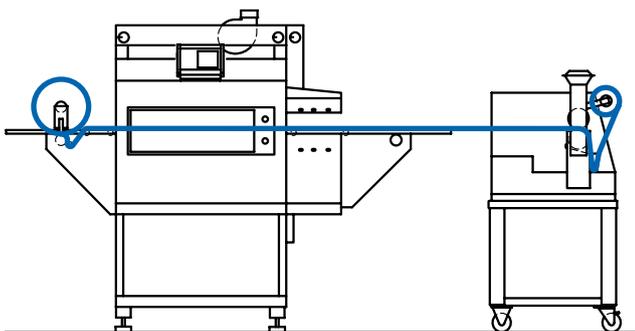

Continuous process flow in combination with MATHIS padder horizontal, vertical type «HVF»

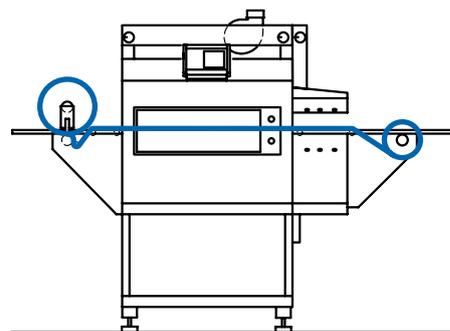

Continuous process flow from roll to roll


HEART OF THE «KTF» - THE NOZZLE FIELD

Discontinuous process flow with specific material holders

Individual material samples are brought up for drying on specific material holders. These then slide through the «KTF» with the prescribed process parameters. The material holders can also slide through the «KTF» reversing – thus the pick-up and deposit point is located at the same place.



The nozzle field, the effective drying area of the «KTF» is divided into an upper and a lower nozzle section. Two regulating flaps can be adjusted so that an even or different ventilation of the upper and the lower nozzle section is achieved. The optimum drying capacity is achieved through a powerful air circulation in combination with a perfectly coordinated heating. The temperature control is based on the resistance measurement. To achieve a precise temperature measurement, the sensor is attached just in front of the air outlet. An additional temperature sensor acquires the temperature of the heating element and protects this against overheating. The heating is dimensioned so that the required drying temperature is reached and maintained quickly.

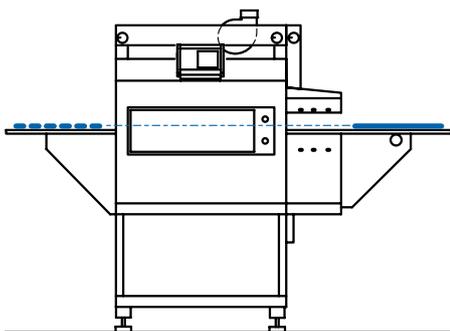

Discontinuous process flow with specific material holders



OPTIONAL, PROCESS-ORIENTATED ADD-ON COMPONENTS

Infrared predryer

The infrared predryer is mounted in front of the dryer field of the «KTF». IR radiators from above and below enable an optimum heat treatment predrying of the material sample. The capacity of the infrared predryer is infinitely variable.

Double field dryer

If different treating temperatures are required (drying, thermolising) or if an increased dryer output must be achieved, several dryer fields are placed one behind the other.

Exhaust fan

Special exhaust fans for the application in heat treatments with powerful air load are available.

Infrared radiation pyrometer

For the contactless measuring of the material sample surface temperature. The IR radiation pyrometer enables the exact and current determination of the surface temperature and from this the resulting residual moisture.

Its emission factor can be adapted individually to the sample to be treated via the «UNIVISION Touch S» process controller.

Material holders in various designs

Material holders such as pin frames adjustable in length and width, pin frames adjustable in length and sieve frames are available for the standardized sample sizes.



Infrared predryer for optimum heat pretreatment



Double field dryer for different treating temperatures, e.g. for drying - thermolising



Easily accessible and simply operated unwinding device



Exhaust air system with adjustable exhaust air zones



Upwinding device, designed as rising roll batcher



«UNIVISION TOUCH S» - THE USER-FRIENDLY PROCESS CONTROLLER

The programmable «UNIVISION Touch S» process controller is optimally adapted to the requirements of its user and to the functionality of the «KTF». It is very simple in operation, thanks to:

- Clearly laid out graphic illustration of the process with all functions.
- Simple and easily remembered function pictograms.
- Menu-driven parameter entry.
- Updating and display of all process data during the operation.
- Information displays (error and fault displays) in plain text.
- Call up of a comprehensive helpfile.
- Update-capable software.
- Import and export of process data to other «UNIVISION Touch S» or to PC.

The «KTF» in combination with the «UNIVISION Touch S» is ready for the present as well as future processes in your business.

OPTIONAL FUNCTIONS OF THE «UNIVISION TOUCH S»

Process programmer

Recurring processes consisting of several treatment intervals can be saved and recalled when required.

Day and week programmer


Not only recurring individual processes, but also complete day and week programs can be predefined and started. For example, the «KTF» can already be preheated at the start of the week before work begins.

Process management using PC

It is possible to connect the «UNIVISION Touch S» with a LAN network using Interface RS232C or Ethernet Interface TCP/IP connection. Process data as well as complete programs are thus created centrally and transferred to the «UNIVISION Touch S».

Even a running process can be changed this way.

Over a connected PC e.g. reports of the processes of the «UNIVISION Touch S» can be created for ISO 9001.

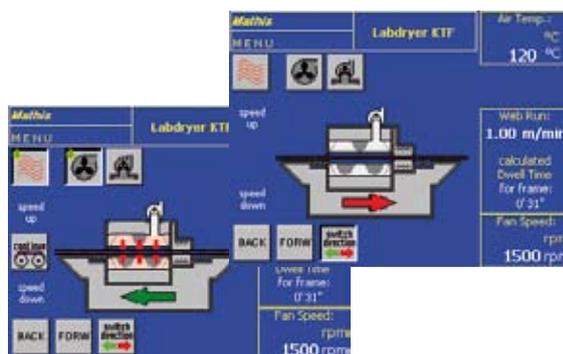


UNIVISION TOUCH S» the clearly laid out, user-friendly process controller

TECHNICAL INFORMATION

	Type «KTF» 350	Type «KTF» 500
Working width:	300 mm	450 mm
Material speed:	0.1 – 2 m/min	
Dwell times:		
Roll – roll	20 – 420 s	
pin frame	15 – 280 s	
Circulation air speed:	up to 6 m/s	
Circulation air capacity:	up to 1000 m ³ /h	
Exhaust air capacity:	up to 400 m ³ /h	
Temperature range:	up to 250 °C	
Heating capacity:	9 kW	13.5 kW
Infrared predryer:	1.8 kW (above the material web) 3.6 kW (above/below the material web)	
Nozzle field length:	800 mm	
Sample sizes:		
width x length	330 x 500 mm	480 x 500 mm
Winding diameter:	max. 250 mm	

«UNIVISION Touch S»

Operating system:	Windows CE
Storage capacity:	16 MB RAM
Touchscreen:	5.7" TFT
Interfaces:	TCP/IP connection (Ethernet Interface), USB, MMC-drive

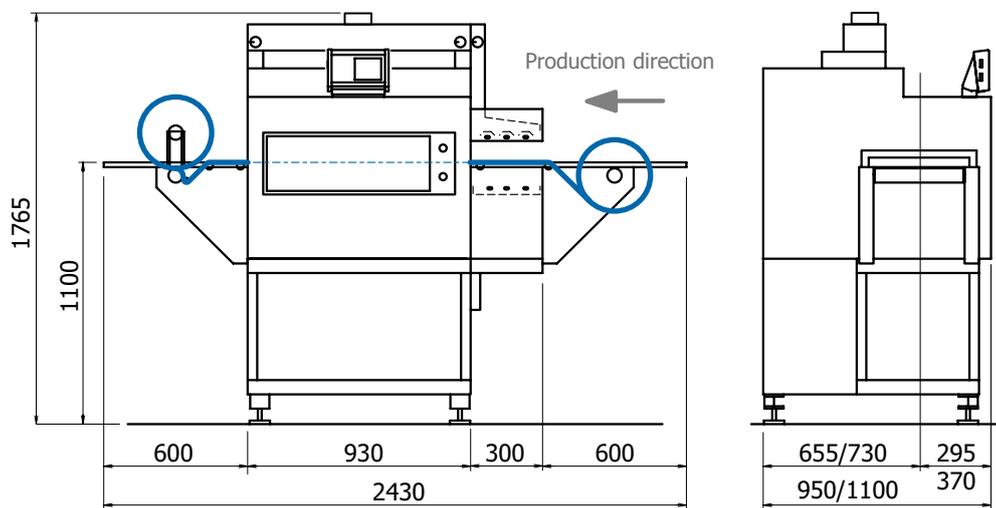

«KTF-S» COATING AND DRYING IN A CONTINUOUS PROCESS

If in addition to the flexible hot air drying, a continuous doctor knife coating is also required, the «KTF-S» comes into use. Naturally, we also plan and optimize the application of the «KTF-S» for you.



OTHER DOCUMENTATION

MATHIS Padder horizontal, vertical type «HVF»

Continuous Coater and Dryer type «KTF-S»

Discontinuous Labdryer type «LTE»

Discontinuous Labcoater type «LTE-S»

Werner Mathis AG

Rütisbergstrasse 3
CH-8156 Oberhasli/Zürich, Switzerland
Telefon: 41(0)44-852 5050
Telefax: 41(0)44-850 6707
E-mail info@mathisag.com
www.mathisag.com

Werner Mathis U.S.A. Inc.

2260 HWY 49 N.E./P.O. Box 1626
Concord N.C. 28026, U.S.A.
Phone: 1-704-786-6157
Fax: 1-704-786-6159
E-mail: usa@mathisag.com

Aparelhos de laboratorio Mathis Ltda

Estrada dos Estudantes 261
CEP 06700 COTIA S.P., Brasil
Phone: 55 (0)11-4612-2333
Fax: 55 (0)11-4612-2598
E-mail: brasil@mathisag.com

Mathis (India) PVT. LTD.

203, Swastic plaza
Pokhran Road No. 2, Thane 400 601
Phone: +91 22 2585 4304
E-mail: india@mathisag.com